





MEGASAN

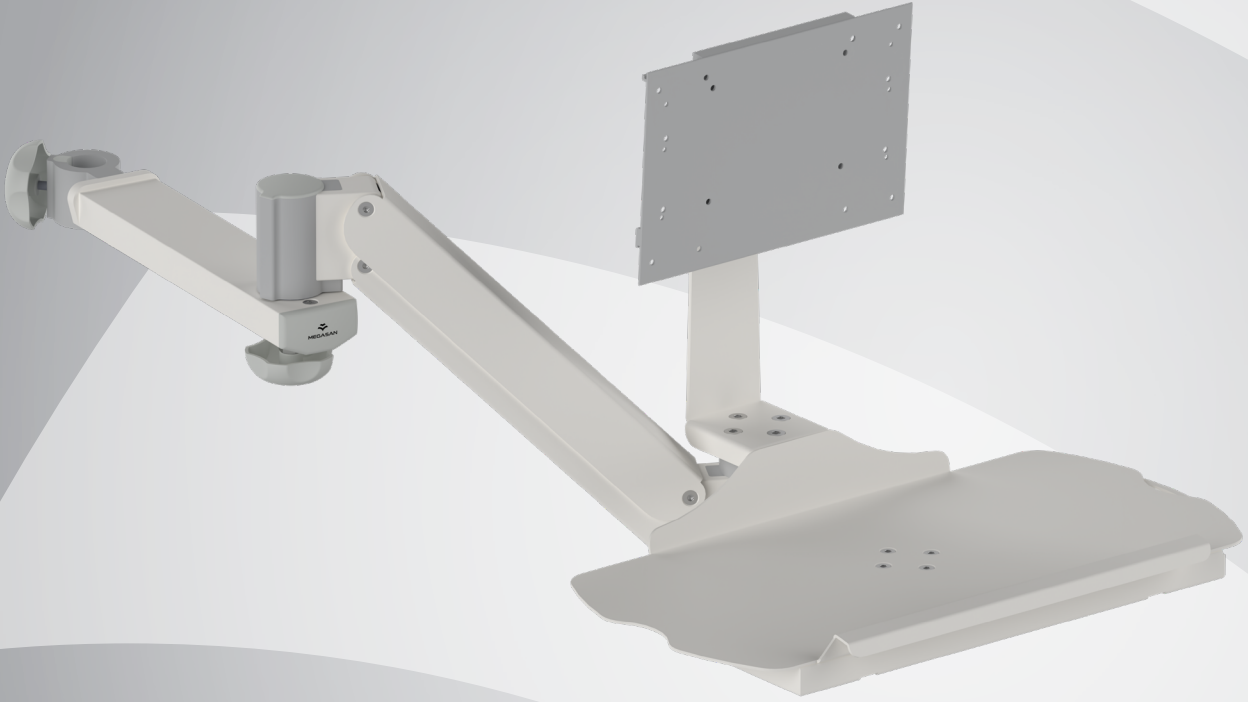

Monitor Stands

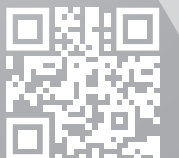

www.megasan.com.tr



MEGASAN Medikal entered the sector in 1981 with assembly and maintenance works and started production in 1998.


It continues its activities in its production facility established on a total area of 70 thousand m² with a closed area of 14 thousand m². There is also a showroom in Ankara where the R&D and sales departments are located. Today, **MEGASAN Medikal A. Ş.** is one of the largest medical gas manufacturers in Europe. In addition to monitor stands in its product range, it produces medical gas plants, area valve service units, medical gas outlets, patient service units, intensive care and emergency service units, pendants, gas adapters.

Today, **MEGASAN** has the rightful pride of offering dozens of best quality products and serves in gas installation projecting, installation and maintenance with its modern manufacturing building equipped with sensitive technological devices and equipments, engineers, extremely well-trained technicians, assembly and maintenance personnels. **MEGASAN Medikal**, is proud of its values and determined mind about getting better day by day.

Our aim; By closely following scientific and technological developments and following a policy parallel to the globalization process in the world, we aim to maintain our leading position in our sector and become a more respected company around the world, while also strengthening this leadership.

Megasan, has the CE (Conformity of Europe) certificate which is very rare in the medical sector in Turkey, proudly states once again that while providing the best service to all its customers;

"The quality of our products is our best reference."

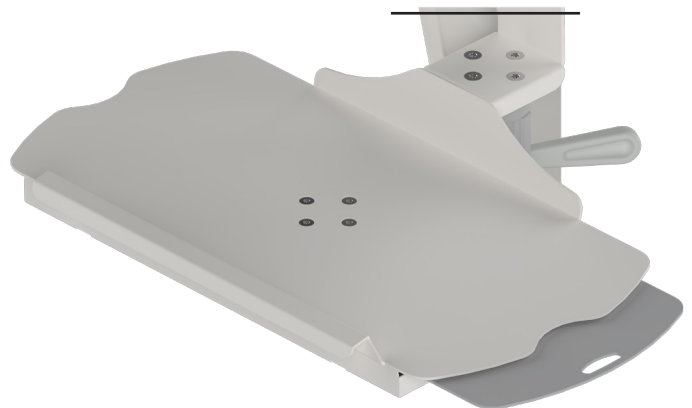



Monitor stands; Patient monitors are equipment with various features that allow medical devices and equipment to be connected to the wall or devices.

Intended Use:

Monitor stands;


The head part has a mechanism that can move in multiple directions,
Allows the monitor to move up and down,
Thanks to its slide mechanism and spring lock system, the position of the monitor can be easily adjusted,


In addition, pendant, intensive care, bedside, etc. where the monitor must be easily attached and removed. It can be used functionally in units.

Mounting holes are drilled in the sled specifically for the monitor.

Environments of Use:


In general, connection equipment can be used in operating rooms, intensive care, service rooms and intervention rooms.

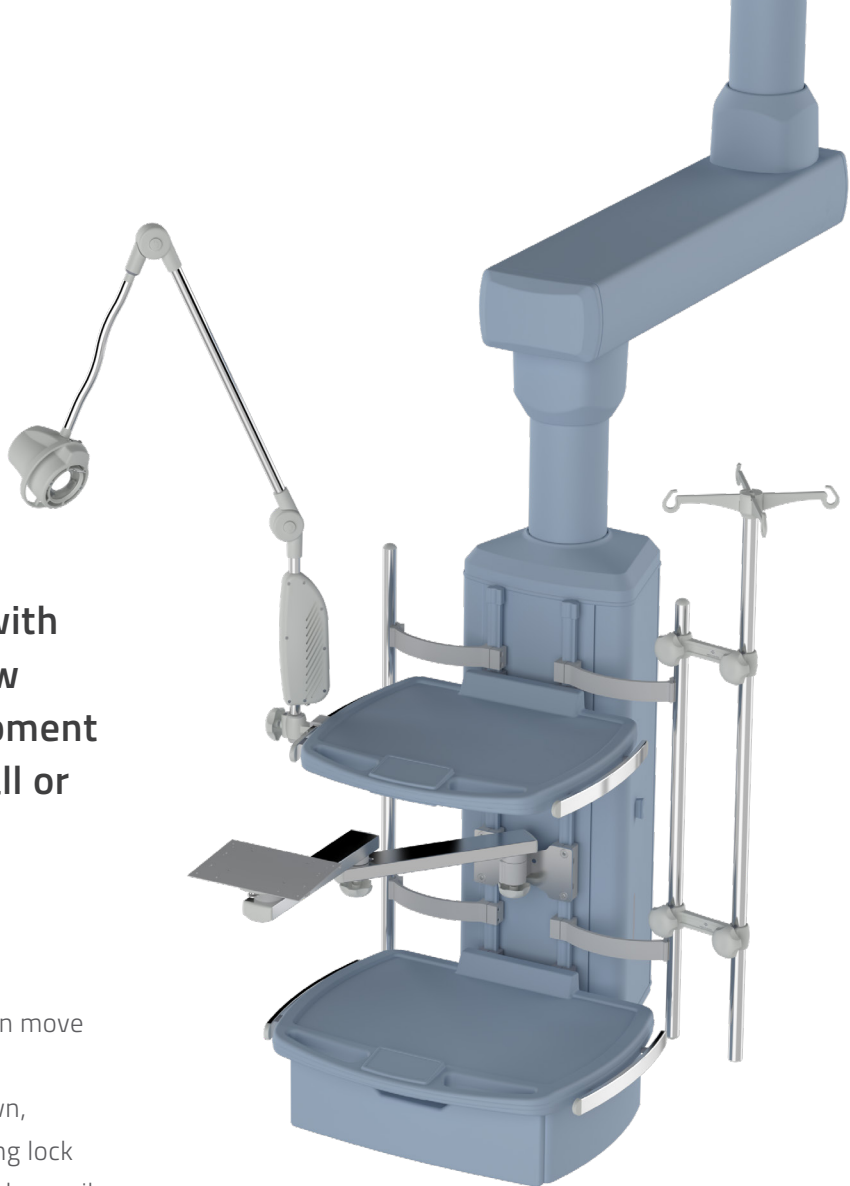

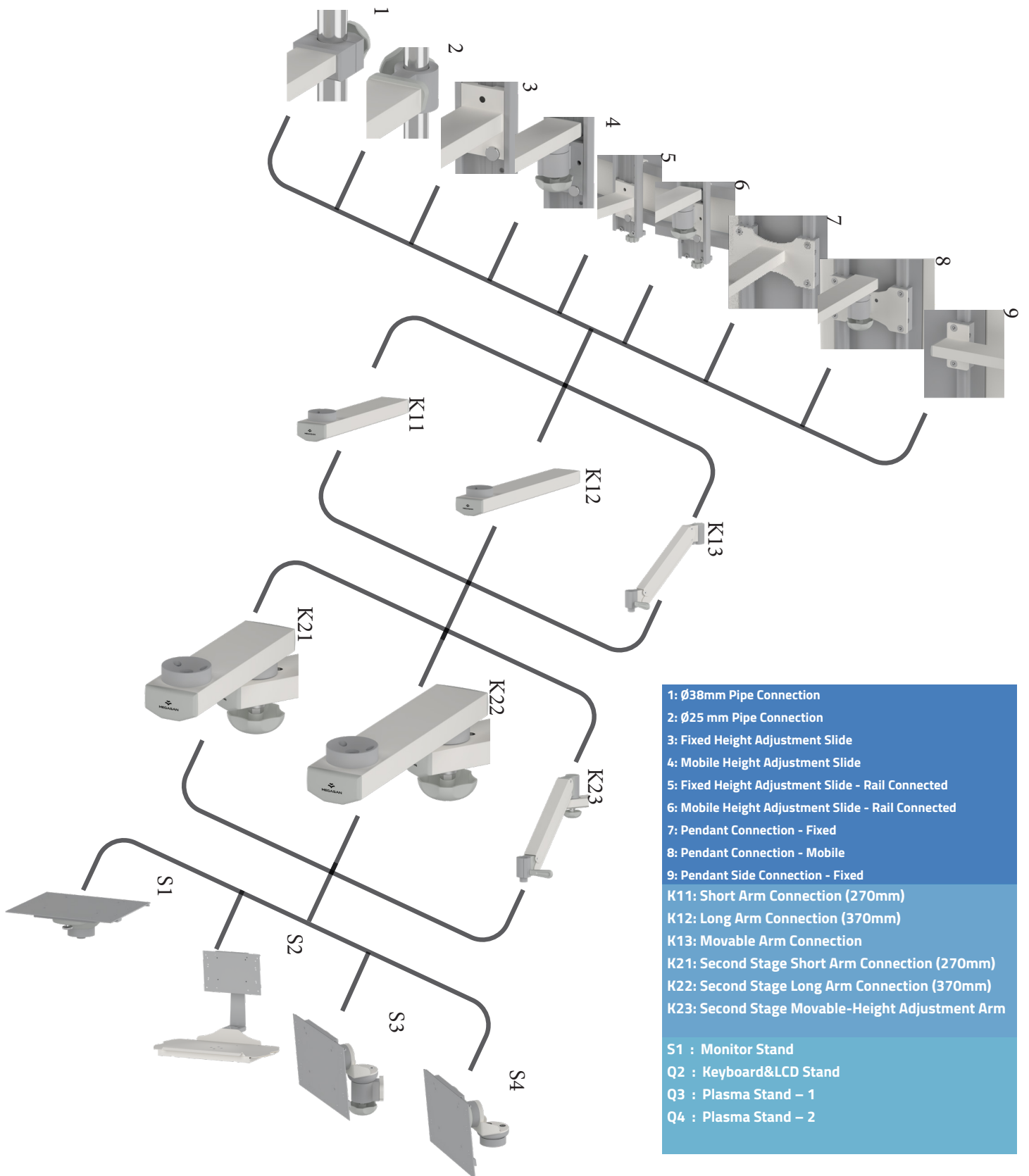

- 1: Ø38mm Pipe Connection
- 2: Ø25 mm Pipe Connection
- 3: Fixed Height Adjustment Slide
- 4: Mobile Height Adjustment Slide
- 5: Fixed Height Adjustment Slide - Rail Connected
- 6: Mobile Height Adjustment Slide - Rail Connected
- 7: Pendant Connection - Fixed
- 8: Pendant Connection - Mobile
- 9: Pendant Side Connection - Fixed
- K11: Short Arm Connection (270mm)
- K12: Long Arm Connection (370mm)
- K13: Movable Arm Connection
- K21: Second Stage Short Arm Connection (270mm)
- K22: Second Stage Long Arm Connection (370mm)
- K23: Second Stage Movable-Height Adjustment Arm
- S1 : Monitor Stand
- Q2 : Keyboard&LCD Stand
- Q3 : Plasma Stand – 1
- Q4 : Plasma Stand – 2

- The table profiles are made of aluminum and anodized for a sleek look.
- Aesthetic and complementary plastic parts provide an elegant appearance to the unit.
- Corrosion-resistant materials are preferred for all parts.
- The monitor connection part can rotate 360 degrees and tilt forward and backward up to 20 degrees.
- Monitors can be easily attached and detached using a sliding profile.
- Monitor tables are available in single and double-arm configurations.
- The arm connection can be fixed or movable as per preference.
- Table heights are fixed, but a movable arm connection should be preferred in places where height adjustment is needed.
- Monitor tables can be attached to intensive care units, pendants, wall-mounted units, and rail-mounted units with various connection components.
- The carrying capacity varies depending on arm preferences.
- Rotating arm connection parts can be tightened after adjustment.
- The height adjustment rails allow vertical movement and can be secured at the desired height.
- Stop bolts are located at the upper and lower parts of the height adjustment rails to prevent the monitor support arms from disengaging from the slide.
- Baskets can be attached to wall rails as a special feature.
- Horizontal or vertical monitor connections are available depending on preference.
- Keyboard & LCD tables feature a sliding system file rack. (S2)

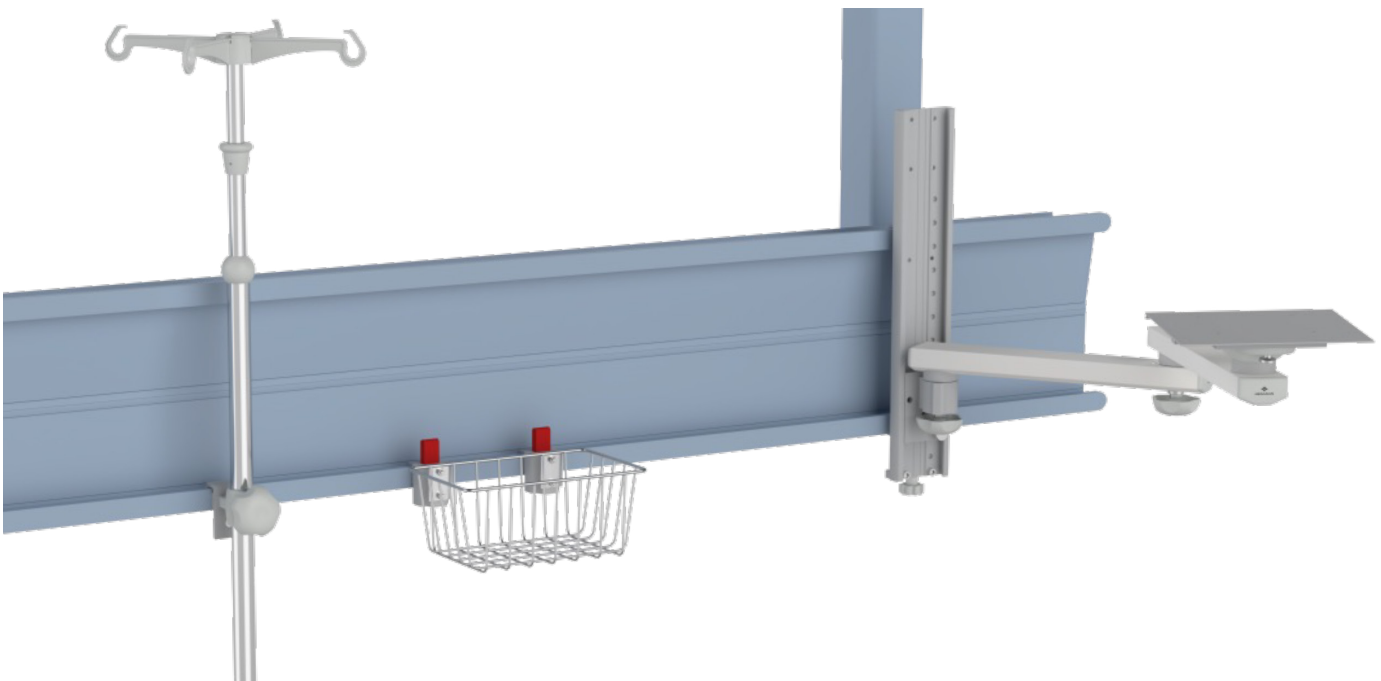


2-K11-S1




2-K12-S1




2-K11-K21-S1

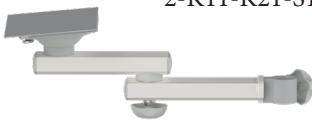

2-K12-K21-S1

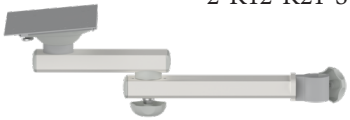

2-K12-K22-S1

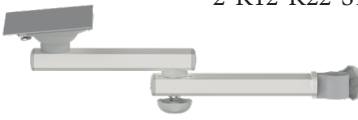


7-K11-S1




7-K12-S1




7-K11-K21-S1

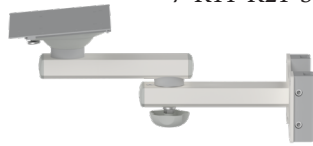

7-K12-K21-S1

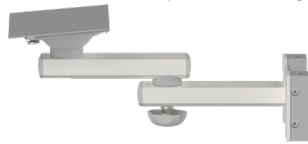

7-K12-K22-S1

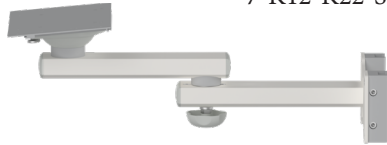


8-K11-S1




8-K12-S1



8-K11-K21-S1

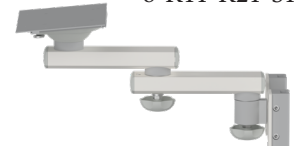

8-K12-K21-S1

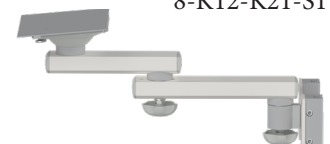

8-K12-K22-S1



5-K11-S1




5-K12-S1




5-K11-K21-S1



5-K12-K21-S1

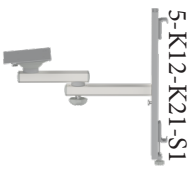

5-K12-K22-S1

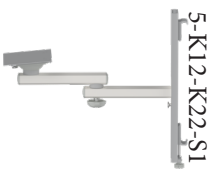

6-K11-S1



6-K12-S1



6-K11-K21-S1



6-K12-K21-S1



6-K12-K22-S1



3-K11-S1-S



3-K12-S1-S



4-K11-S1-S



4-K12-S1-S






2-K11-K23-S2



2-K12-K23-S2

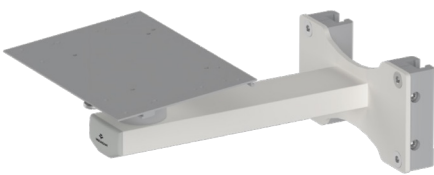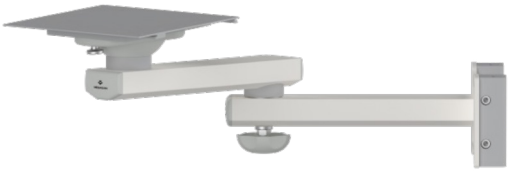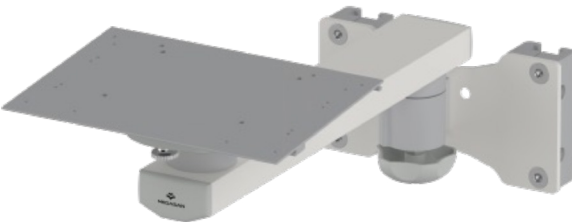

Monitor Stand with Pipe Connection (Toblerone) - Ø38mm	
	
Monitor Stand with Pipe Connection (Single Arm) - Ø25mm	Monitor Stand with Pipe Connection (Double Arm) - Ø25mm
	
Monitor Stand Wall Mounted Basket (Single Arm) Fixed	Monitor Stand Wall Mounted Basket (Single Arm) Rotary
	
Monitor Stand Wall Mounted Basket (Double Arm) Fixed	Monitor Stand Wall Mounted Basket (Double Arm) Rotary
	
Monitor Stand Pendant Mounted Single Arm Fixed	Monitor Stand Pendant Mounted Double Arm Fixed
	
Monitor Stand Pendant Mounted Single Arm Rotary	Monitor Stand Pendant Mounted Double Arm Rotary
	

HARNESSING THE POWER OF THE SUN FOR A SUSTAINABLE FUTURE!



MEDİKAL GAZ SİSTEMLERİ TIBBİ ALET ve CİHAZLAR A.Ş.

SHOWROOM

Yenibahar Mah. 24/12/19 06370 Batıkent Ankara-TURKEY

+90 (312) 865 18 18 - 865 18 21 +90 (312) 865 18 22

HQ

Cumhuriyet Mah. Zafer İş Hanı No:33/411 18100 Merkez/Çankırı-TÜRKİYE

FACTORY

Dere Mahallesi Sağlık Sokak No:29 18700 Eldivan/ÇANKIRI-Türkiye

+90 (376) 218 10 25 - 218 14 07 +90 (376) 218 12 26

+90 (541) 218 10 25 +90 (533) 323 61 63 +90 (555) 982 48 98

DISTRIBUTOR

info@megasan.com.tr www.megasan.com.tr

Rev.2.0-09.10.2023